User Guide

PAP Link

Download on the App Store

GET IT ON Google Play

Install PAP Link

- **For iOS users:**

  Search "PAP Link" in App Store


- **For Android users:**

  Search "PAP Link" in Google Play

Login and register


Login

Using a verrification code or password to login.

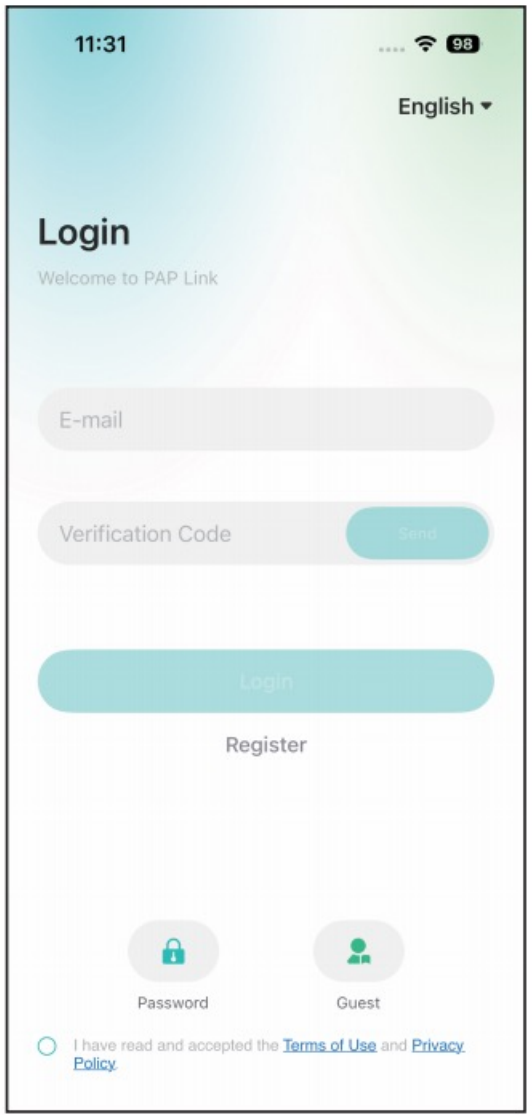

Register

Click "Register". Follow the prompts to complete the registration.

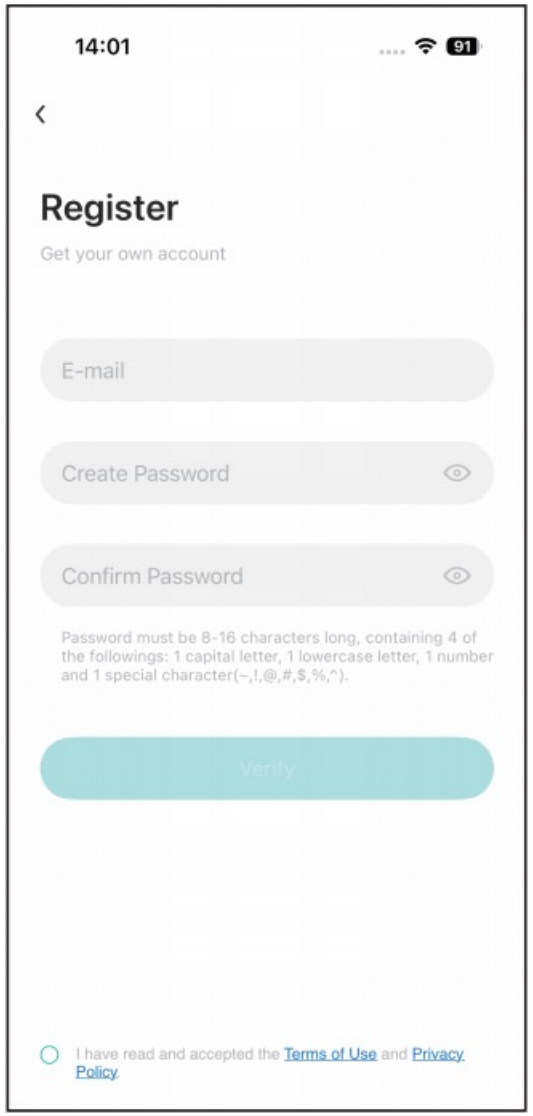


Note:


Read and tick "I have read and accepted the Terms of use and Privacy policy.", otherwise login or registration will fail.

Add Device


Account and Password Login

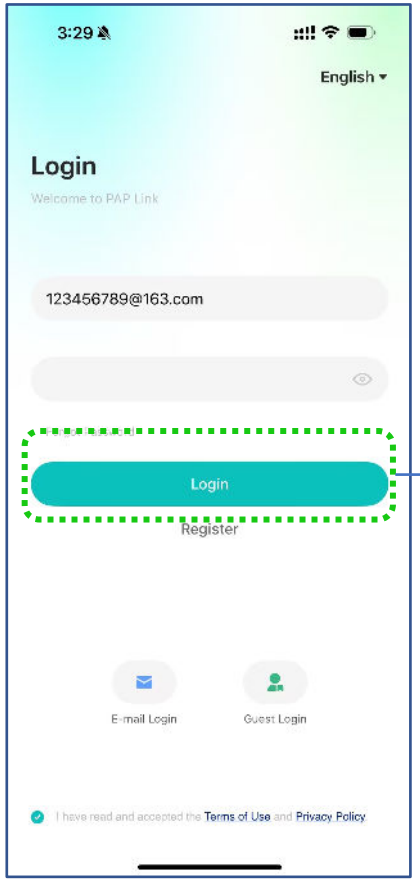


Add Device

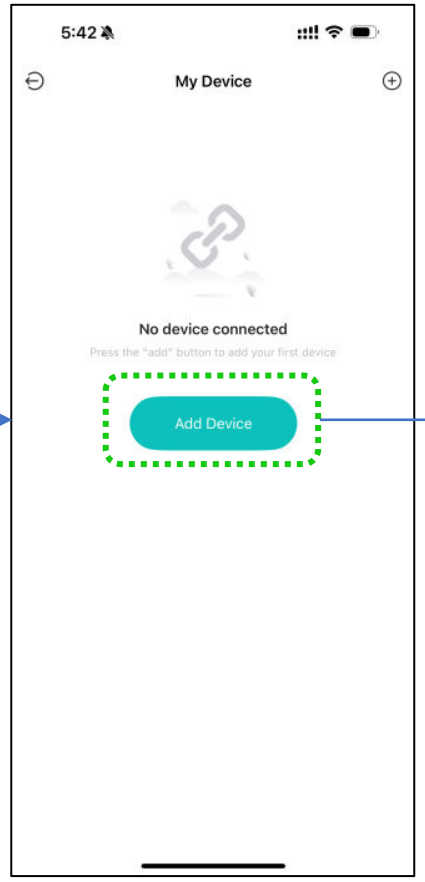

Search and Connect Device

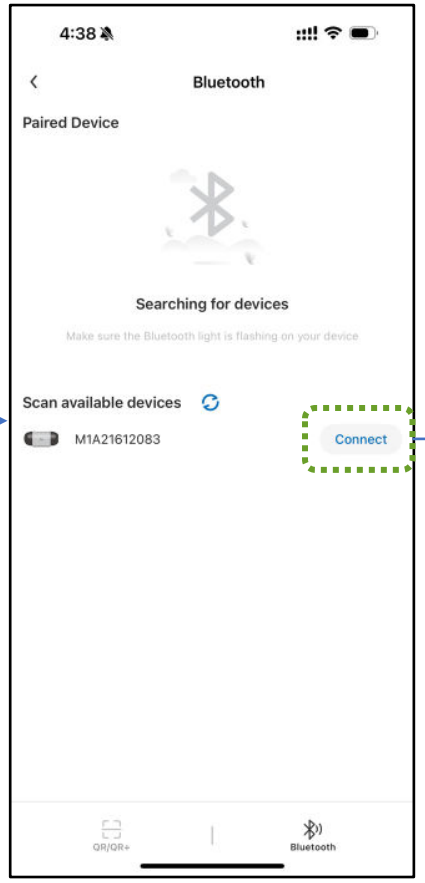


Entered Homepage Successfully

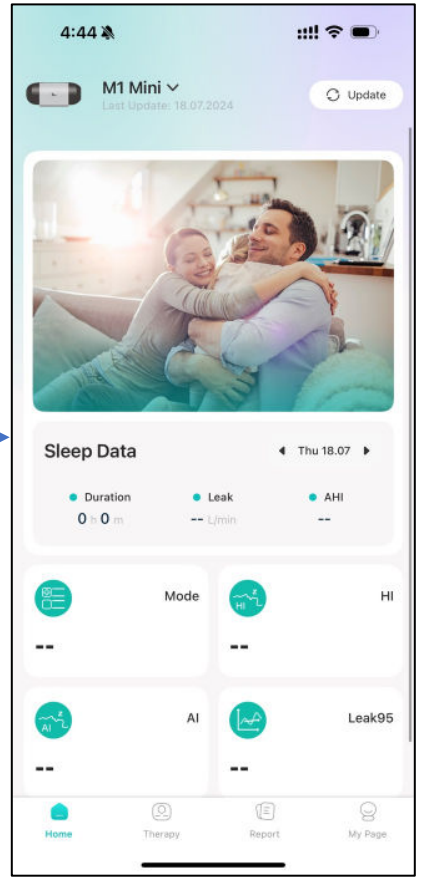


Bluetooth-enabled Models

Add Device


Account and Password Login

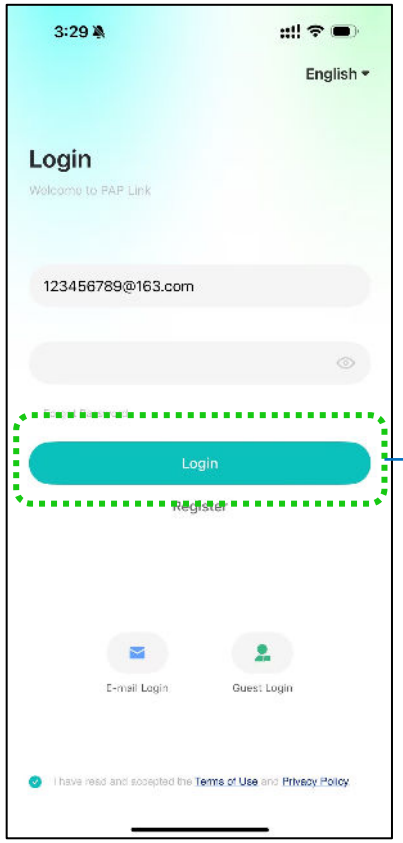


Add Device

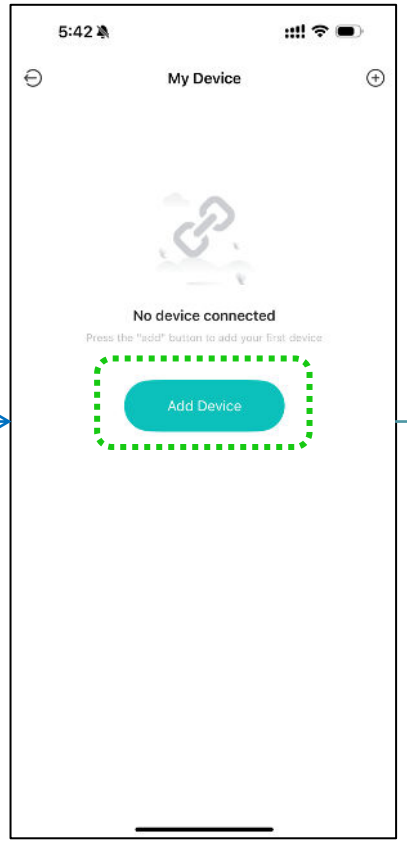


Enter QR/QR+ Scan Interface

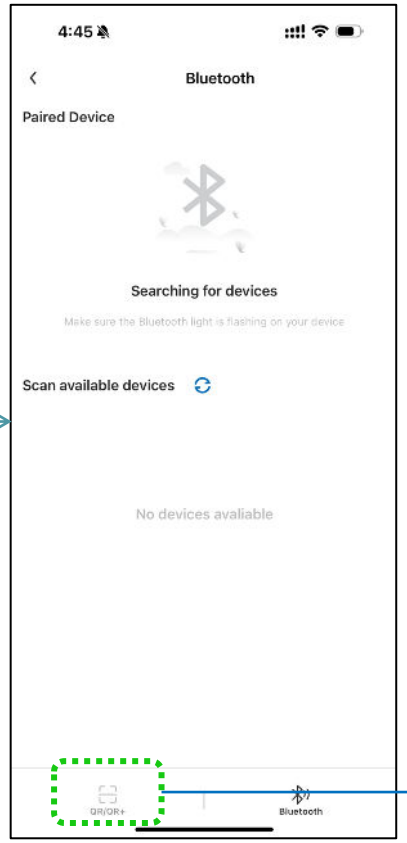

QR/QR+ Scan Interface

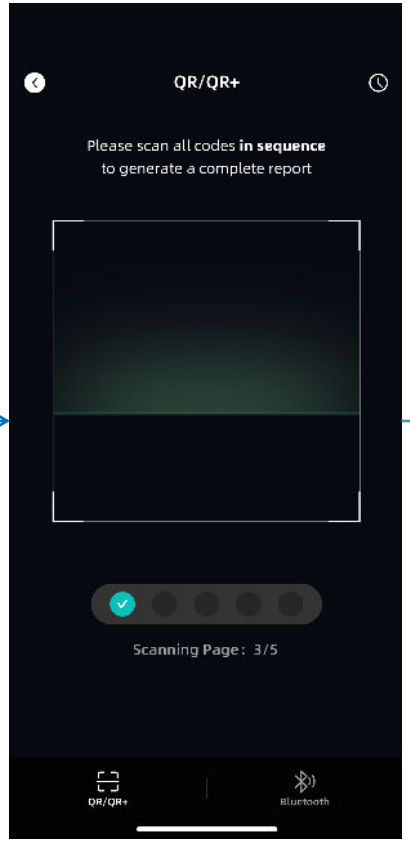


Entered Homepage Successfully

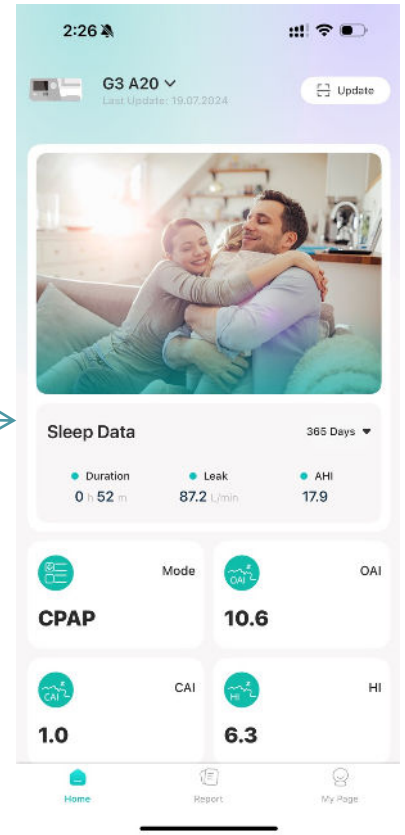


Non-Bluetooth Models - QR/QR+ Scan

Add Device Notes

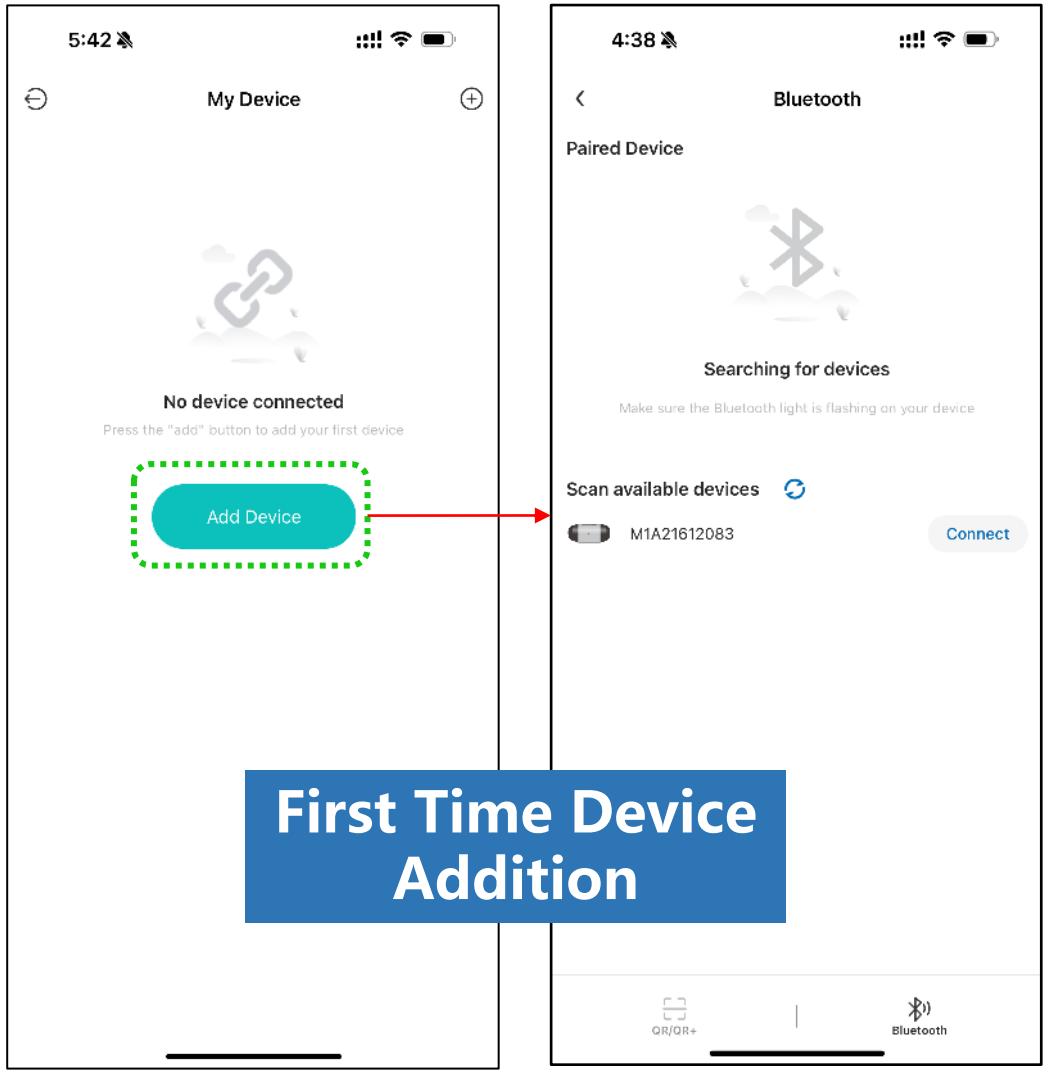

First Time Device Addition

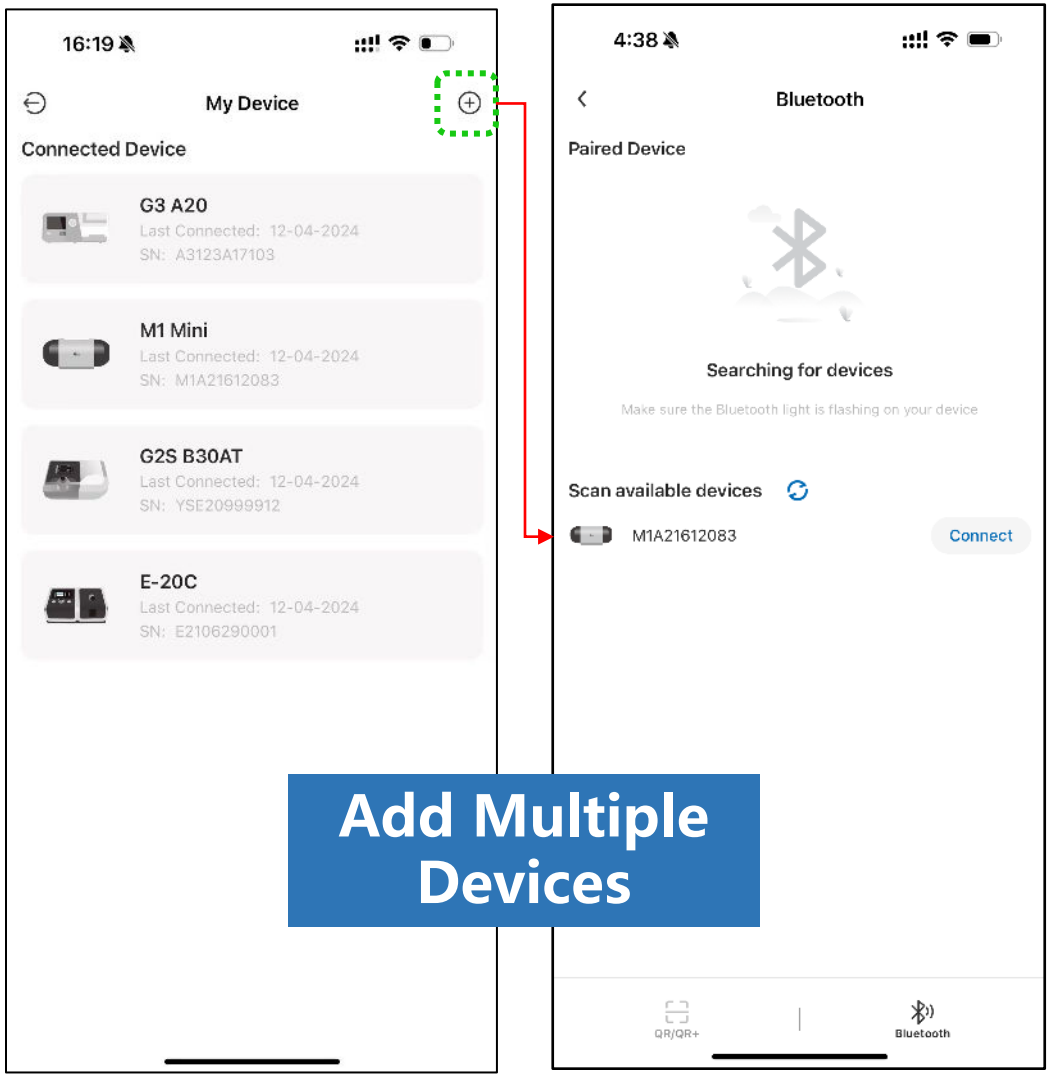



Add Multiple Devices

Add Device Notes

Old Data Structure Devices


Multiple QR+ codes on the device interface cannot be manually switched; the flashlight is used to switch the QR+ code.

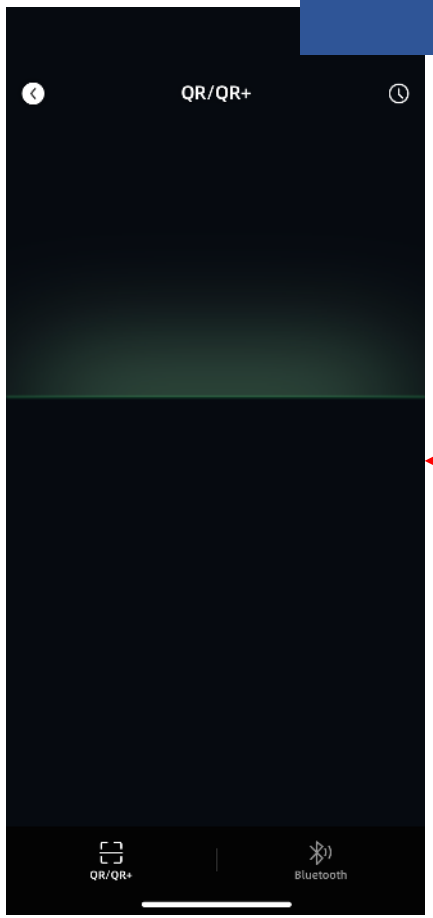

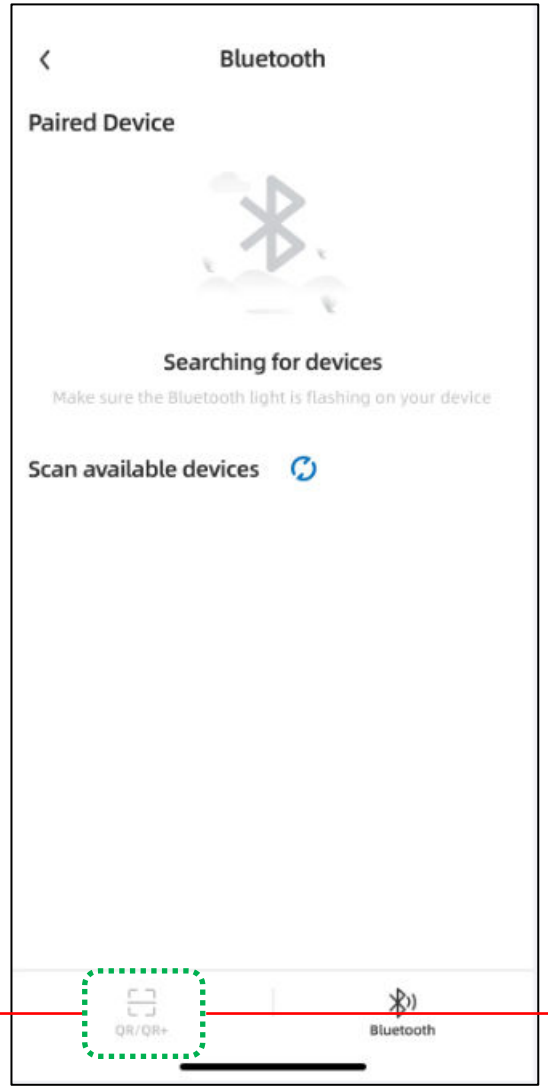

New Data Structure Devices


Multiple QR+ codes on the device interface have been changed to manual switching.

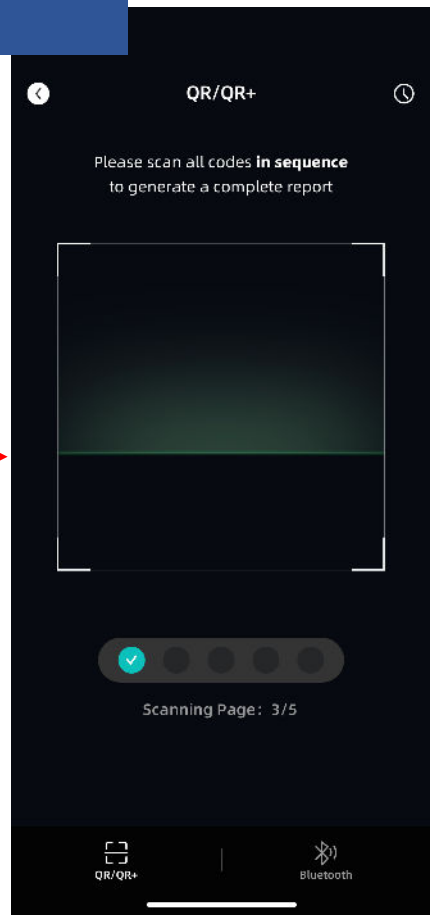

Home interface




Click on the top left corner to access the My Device interface. On My Device interface, the user can add device, delete device and log out.

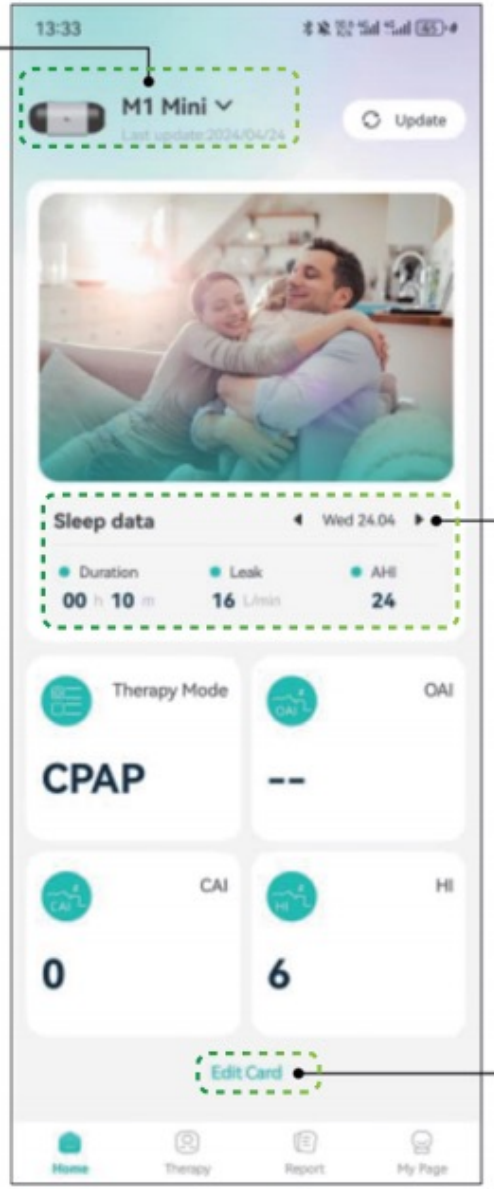

Get Sleep data for a specified date.

Select date

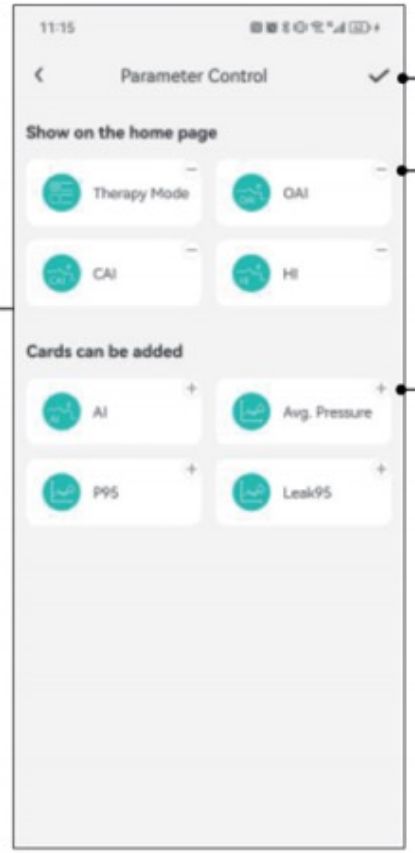


Confirm

Remove


Add

Click “Edit Card” to modify the parameters displayed on the Home interface.


Home interface Notes

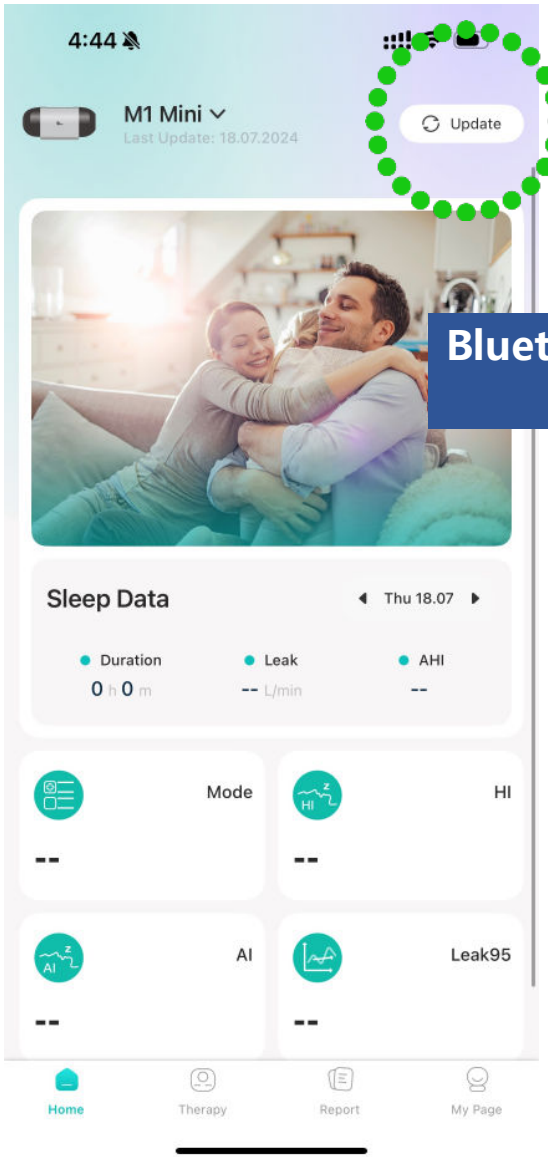


Bluetooth-enabled Models

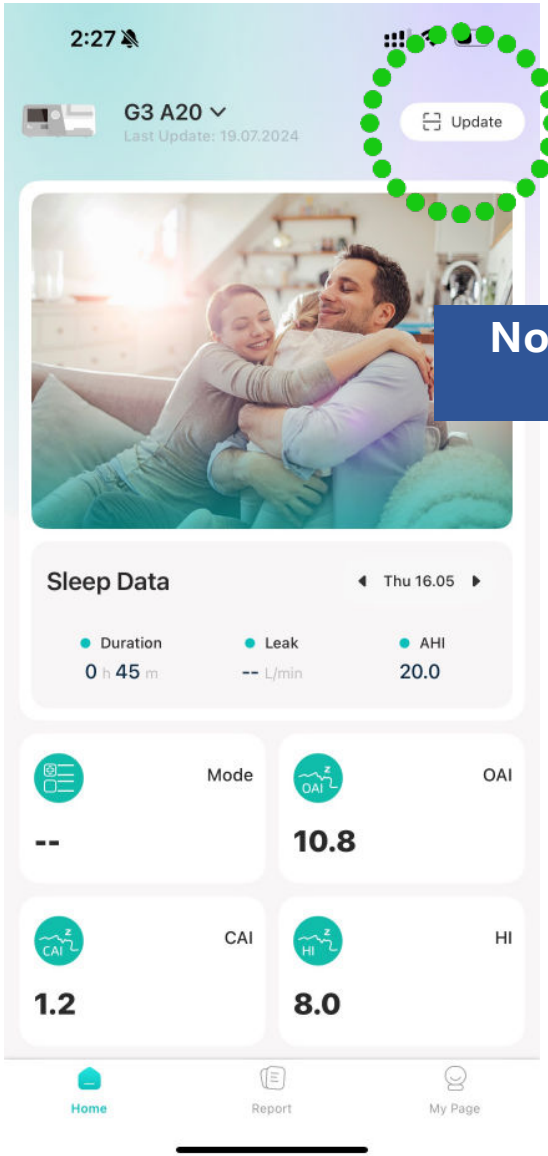


Non-Bluetooth Models

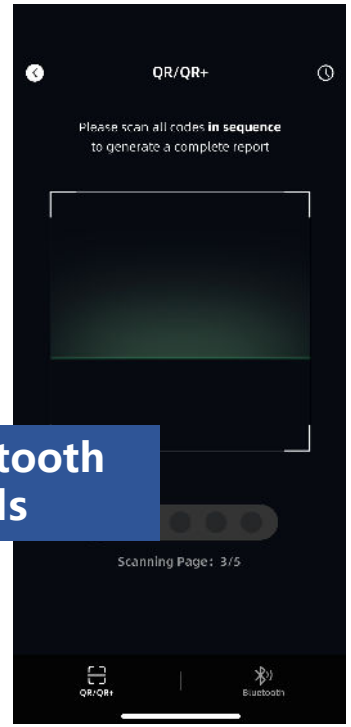

Home interface Notes

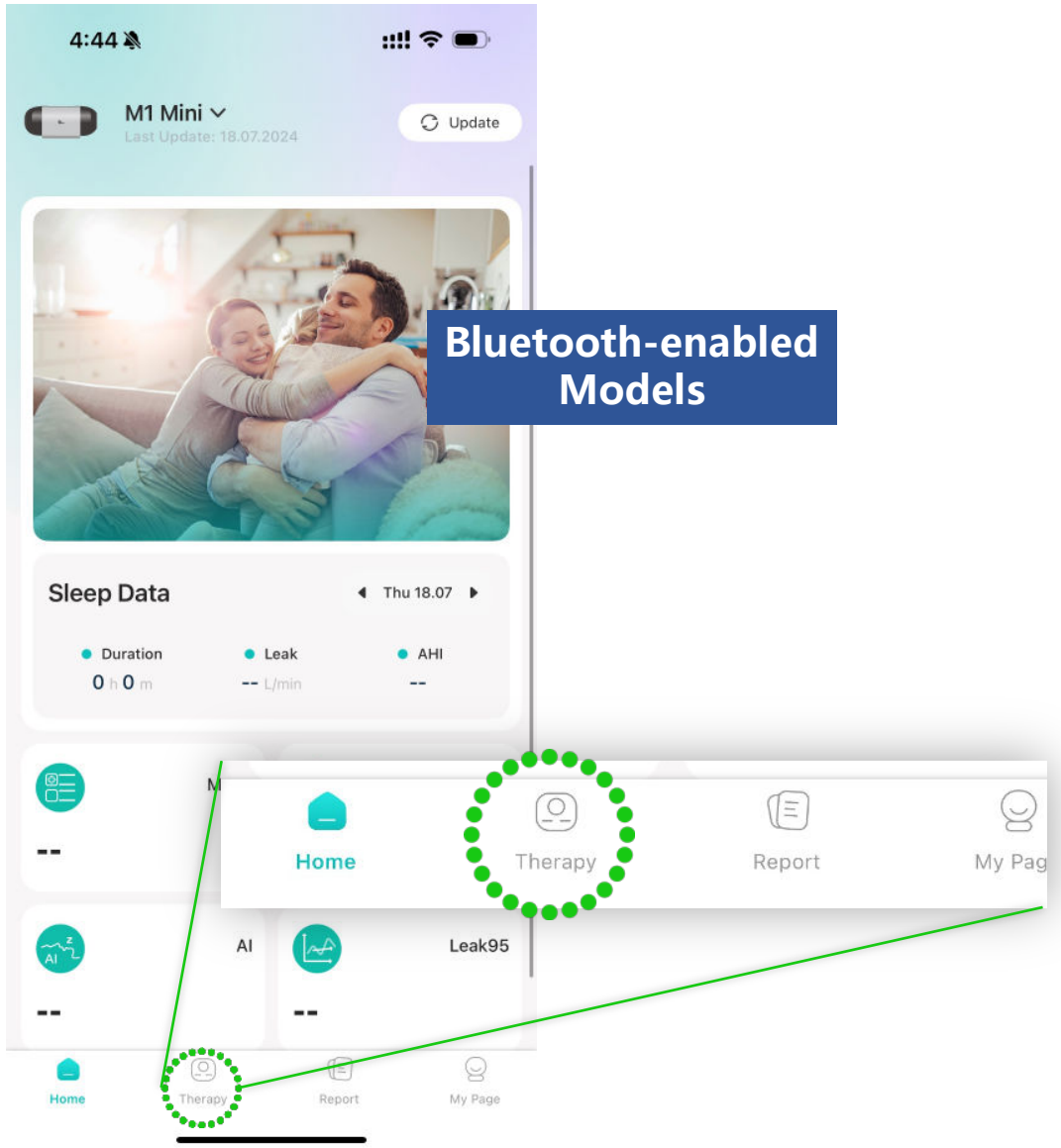


Bluetooth-enabled Models

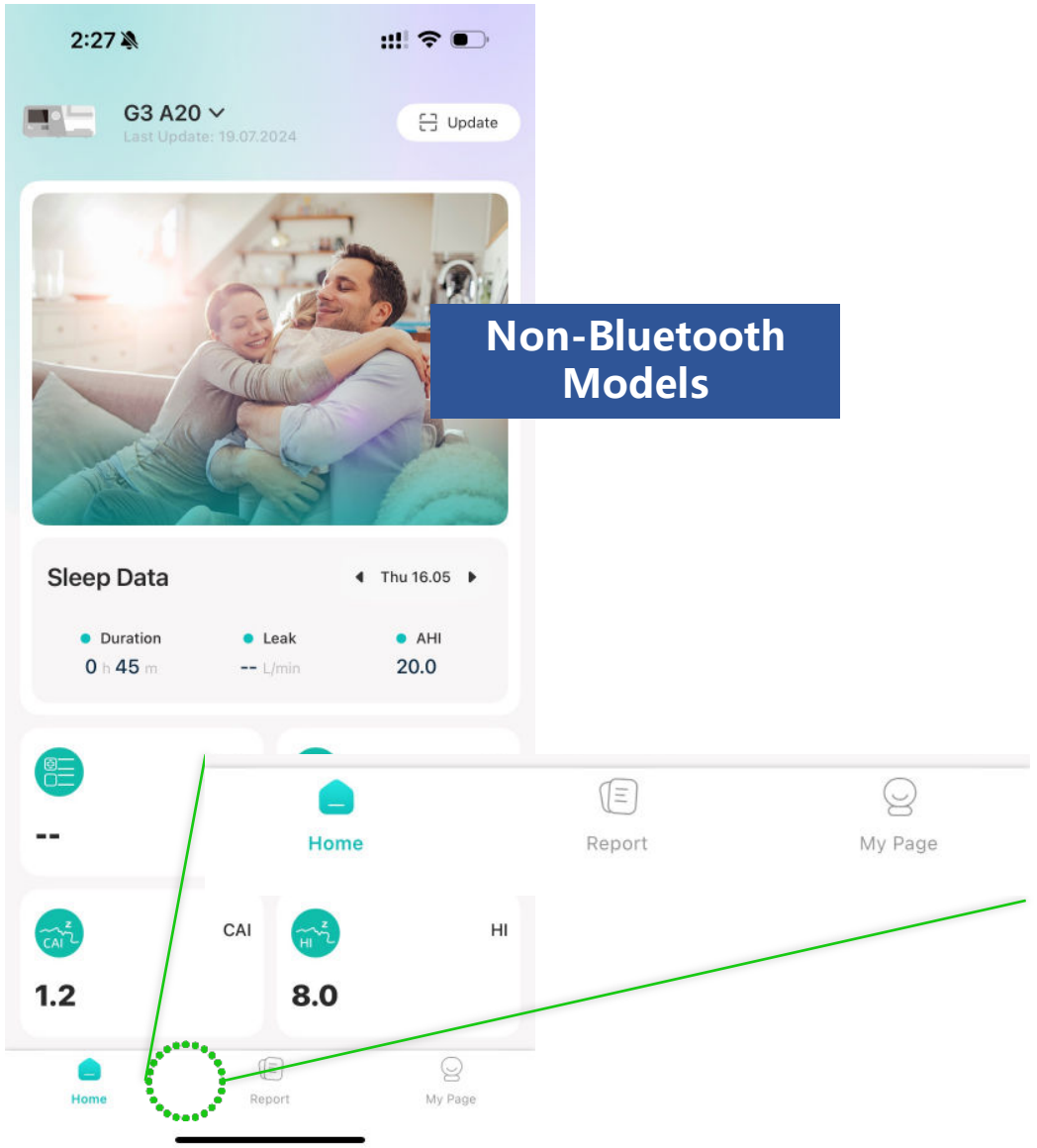

Non-Bluetooth Models


Start Treatment

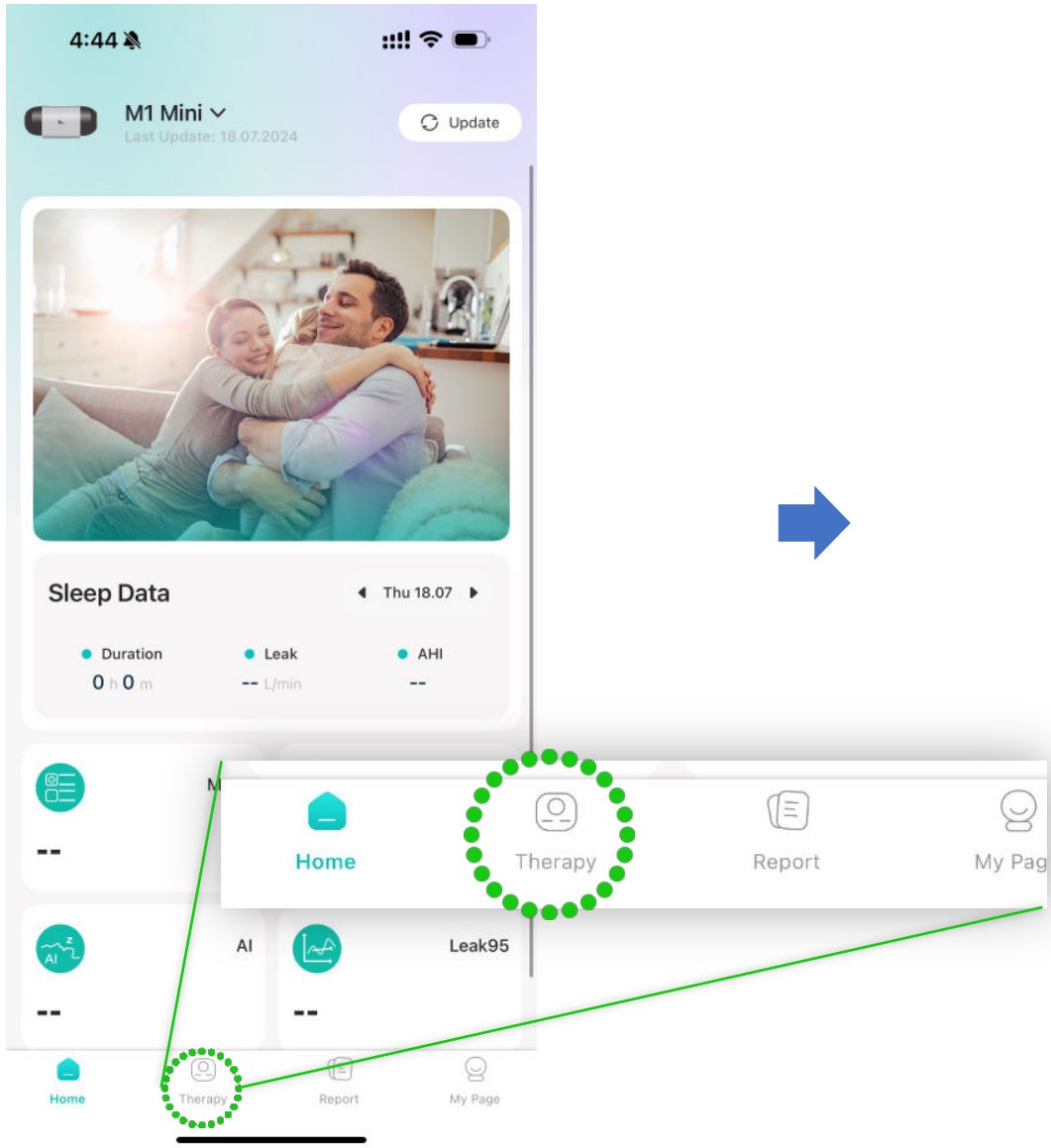

Standby Interface




Post-Start Interface

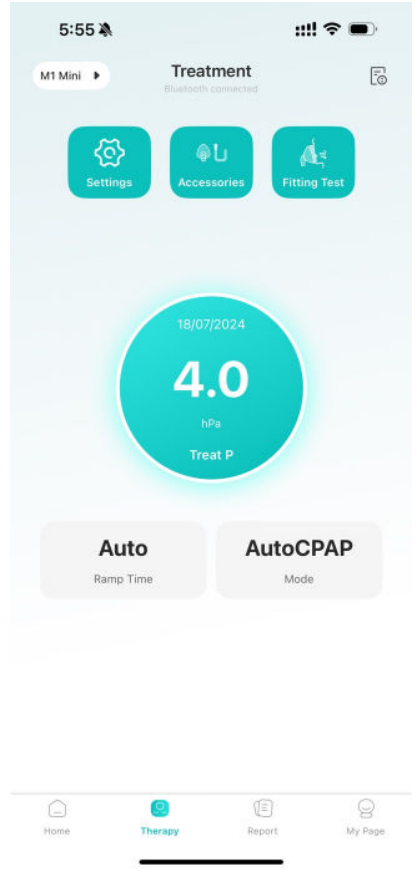


Start Therapy (Bluetooth-Controlled Models Only)

Start Treatment


Enter Patient Settings

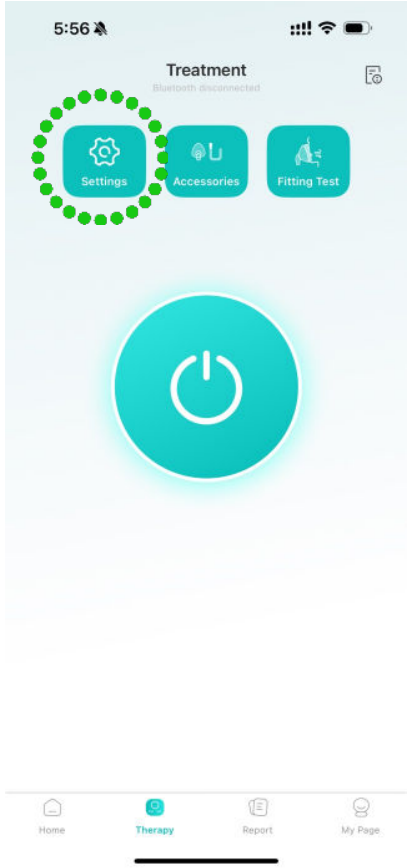


Patient Settings Interface

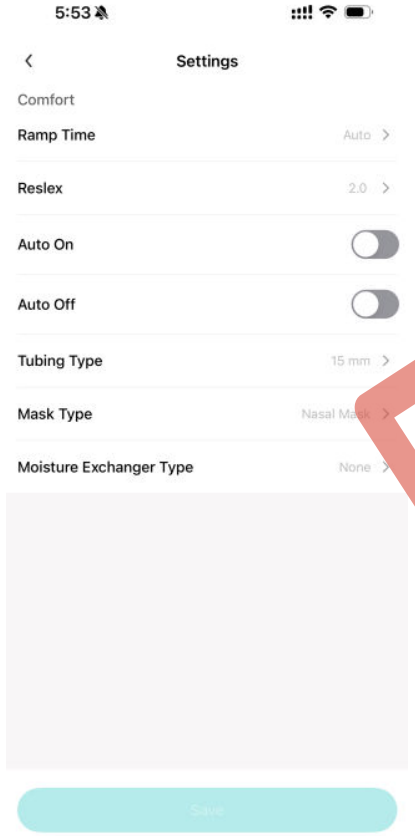

Enter Doctor Settings

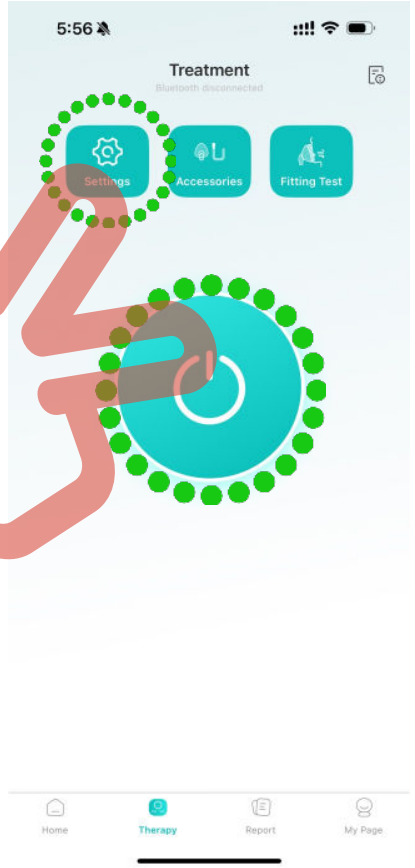

Doctor Settings Interface

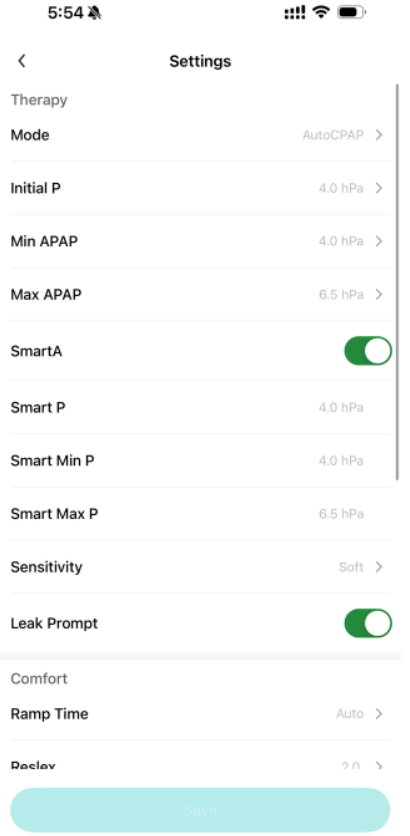


Patient Parameter Adjustment


Doctor's Parameter Adjustment

Get or update report

For the device connected via Bluetooth wireless technology:

1. Click Report to access the Report interface to get report.
2. Upon connection to the PAP Link App, the system will automatically check for any available updates. Alternatively, users can manually trigger the update process by clicking ↻.

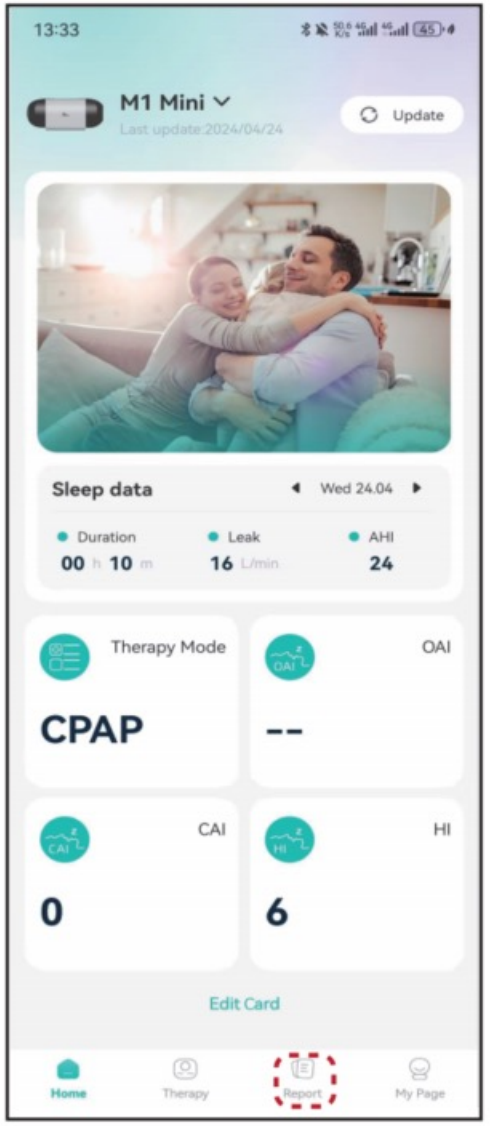

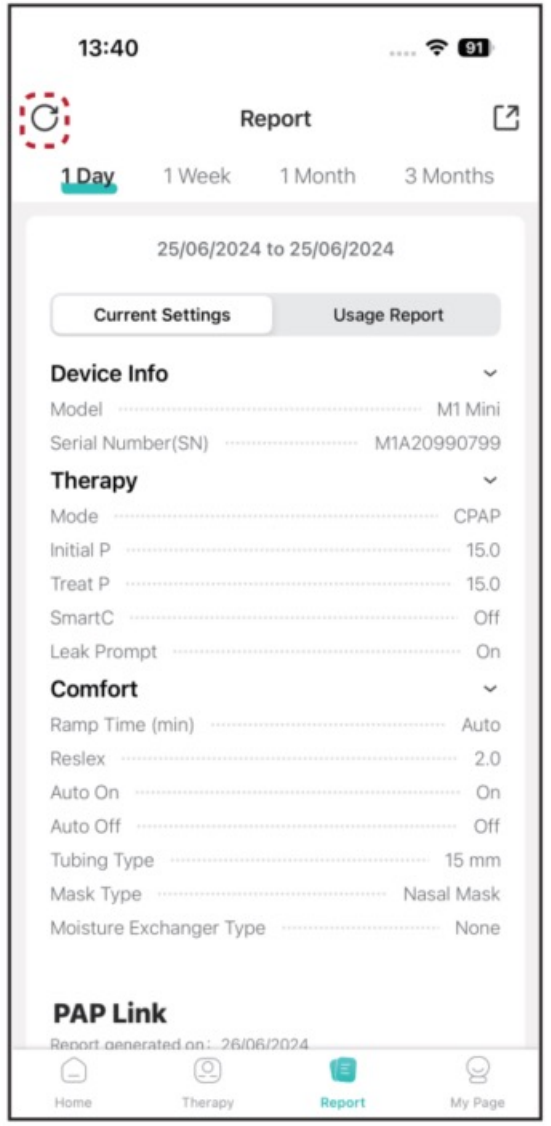

Get or update report

For the device connected via iCode QR/QR+ code scanning:

1. Find the iCode QR/QR+ code on the device.
2. Click Report . Click on the top left corner to scan the iCode QR/QR+ code to get report.
3. Initate a scan to update the report.

iCode QR/QR+





Note: Please note that viewing the report requires an active mobile network connection.

The inmages provided are for reference purposes only and may not accurately depict the actual report interface, which can vary with the specific device model.

Share report

Click  to share the report.

The user can choose the preferred sharing method to distribute the report to other applications.

